

SE 927

Installation, Maintenance and Operating Instructions

DESCRIPTION

Bench mounted AERATED LABORATORY EYE WASH unit, squeeze handle stay open activation valve with dual aerated eye wash heads. Supplied with Emergency Eye Wash sign. **Manufactured to comply with AS 4775-2007.**

SPECIFICATIONS

EYE WASH	Aerated, Dual Stream
VALVE	10mm Chrome Plated Brass
ACTIVATOR	Squeeze Handle

INSTALLATION

After unpacking carton, ensure all components are present. Mark out the desired position on the bench top (that has a maximum thickness of 57mm). Drill a 25mm - 30mm hole. Mount the unit to the bench top using the hardware supplied. Connect water supply to the unit. Turn water supply on and assure all connections are water tight. Activate the unit and allow to flow for 30 seconds to flush line of any impurities.


SPECIFICATIONS

The Australian AS4775 and ANSI 358.1) states that all safety equipment shall be activated on a weekly basis to flush the line and verify proper operation. All units furnishes a testing tag with each unit. On this tag the date of the inspection and the inspectors initials should be noted.

NOTE : This unit should be connected to an uninterrupted source of potable water, no less than 15mm (1/2"), with a minimum flowing pressure of 130KPa and a maximum static pressure of 860 KPa.

OPERATION

The eye wash is activated by squeezing the activation handle. The unit will operate until the activation handle pin is lifted to turn off the unit.

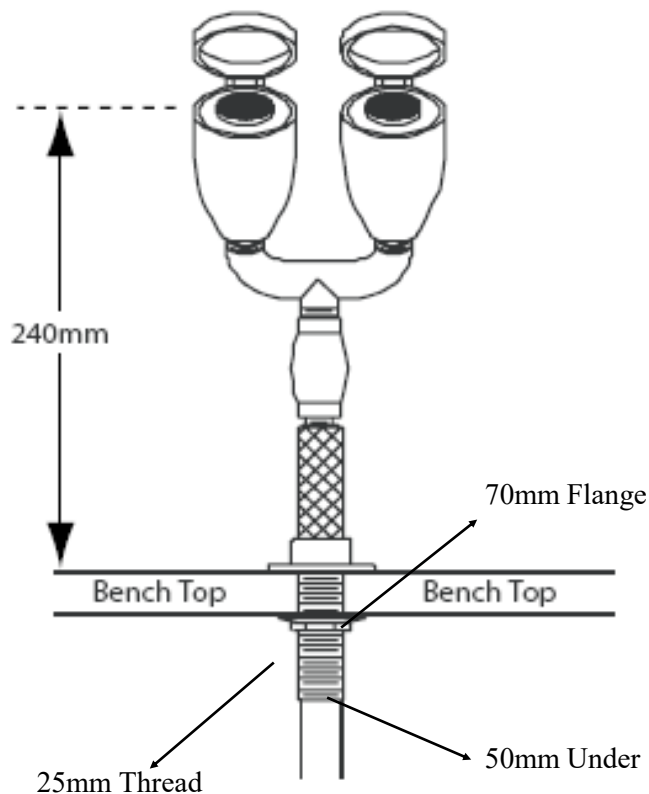

FLOW DATA

FLOW RATE: (LPM): 12

WARNING : Use only genuine Speakman parts when replacing or repairing components. When ordering quote correct part numbers listed.

AUSTRALIAN  SCIENTIFIC

Australian Scientific Pty Ltd
11 McDougall Street
Kotara NSW 2289

E: sales@austscientific.com.au
W: www.austscientific.com.au
Free call: 1800-021 083



SPEAKMAN®
SAFETY PRODUCTS

SE 927
Bench Mount Aerated
Eye Wash

Installation, Maintenance and Operating Instructions

MAINTENANCE

Maintenance should be carried out weekly and yearly as detailed in the schedule. Should you need to repair this unit, please reference the spare parts list for correct part numbers. Before any maintenance is **done be sure the water supply is shut off.**

MAINTENANCE SCHEDULE

WEEKLY	Activate eyewash unit using push handle Allow water to flow for 30 seconds to flush line. Ensure water flows evenly from both outlets. Ensure water is clean.
YEARLY	Ensure unit complies with AS 4775 –2007

1 RPG05 1177—Sold in a packet of 6



3 Complete Eye Wash Assembly Head



8 Optional - Wall Bracket Not supplied With the SE 927




Eye Wash Sign 125 x 245mm supplied



SPARE PARTS

No	Description	Part No.
1	Aerator Screen & Cap	RPG05 1177
2	Yellow Flip Cap Dust Cover	07 0104
3	Aerated Nozzle Head Ass	G38 0415
4	Bull Horn Tee & Nozzles	G24 0110 parts 3 & 4
5	Squeeze Handle & Valve	RPG20 1742
6	Deck Flange, Nut	SE DECKFL
7	1.8m Hose Assembly	G63 0064
8	Wall Mount Bracket S/Steel	SE 927BRK

WARNING : Use only genuine Speakman parts when replacing or repairing components. When ordering quote correct part numbers listed.

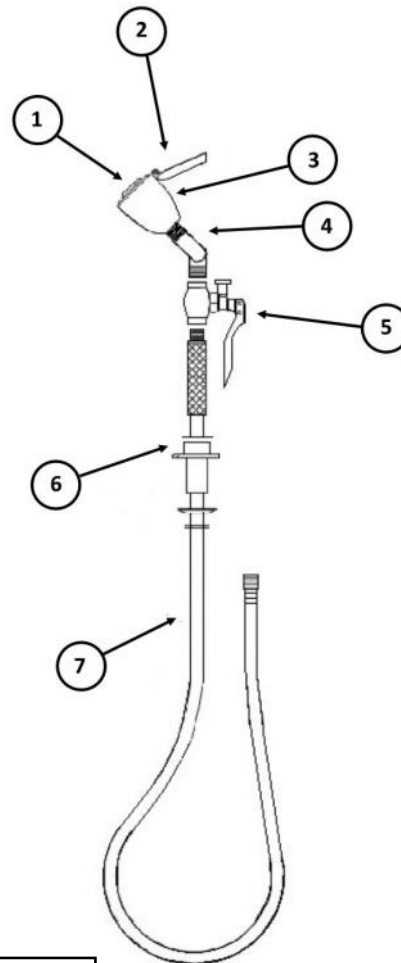

SPEAKMAN®
SAFETY PRODUCTS

SE 927
Bench Mount Aerated
Eye Wash