

AU25302 Poly CW

150 Litre Corrosive Substance Storage Cabinet

With Dividing Wall

for Acids & Bases



AU25302 Poly CW

Polypropylene Corrosive Cabinet

With Dividing Wall

Complies with AS3780

150 Litre Maximum Capacity

External Size: 1000mm H x 1000mm W x 500mm D

Supplied with two shelves as standard.

The cabinets are fully compliant with the relevant Australian Standards.

Fire Retardant Polypropylene construction Corrosive Cabinets stops the effects of rusting when using highly corrosive vapour acids tested to AS 1530.3—1989.

Locking Mechanism, complies to the AS 3780 Australian Standard.

No Internal Metal Parts.

Adjustable and removable shelving within the cabinet.

Each cabinet is provided standard with a liquid tight sealed sump to contain spills or leaks within the cabinet.

Venting position options are included with the cabinet.

Closing of cabinet - spring loaded door hinge with manual slide door handle closing locking.

Opening of cabinet - manual slide handle opening with self closing doors.

Individual key lockable doors.

12 Months Manufacturers Warranty on parts and labour.


Clean with warm soapy water.

Clean up any spills within the cabinet immediately.

After sales service and training available

AU25302 Poly CW

150 Litre Corrosive Substance Storage Cabinet
With Dividing Wall
for Acids & Bases

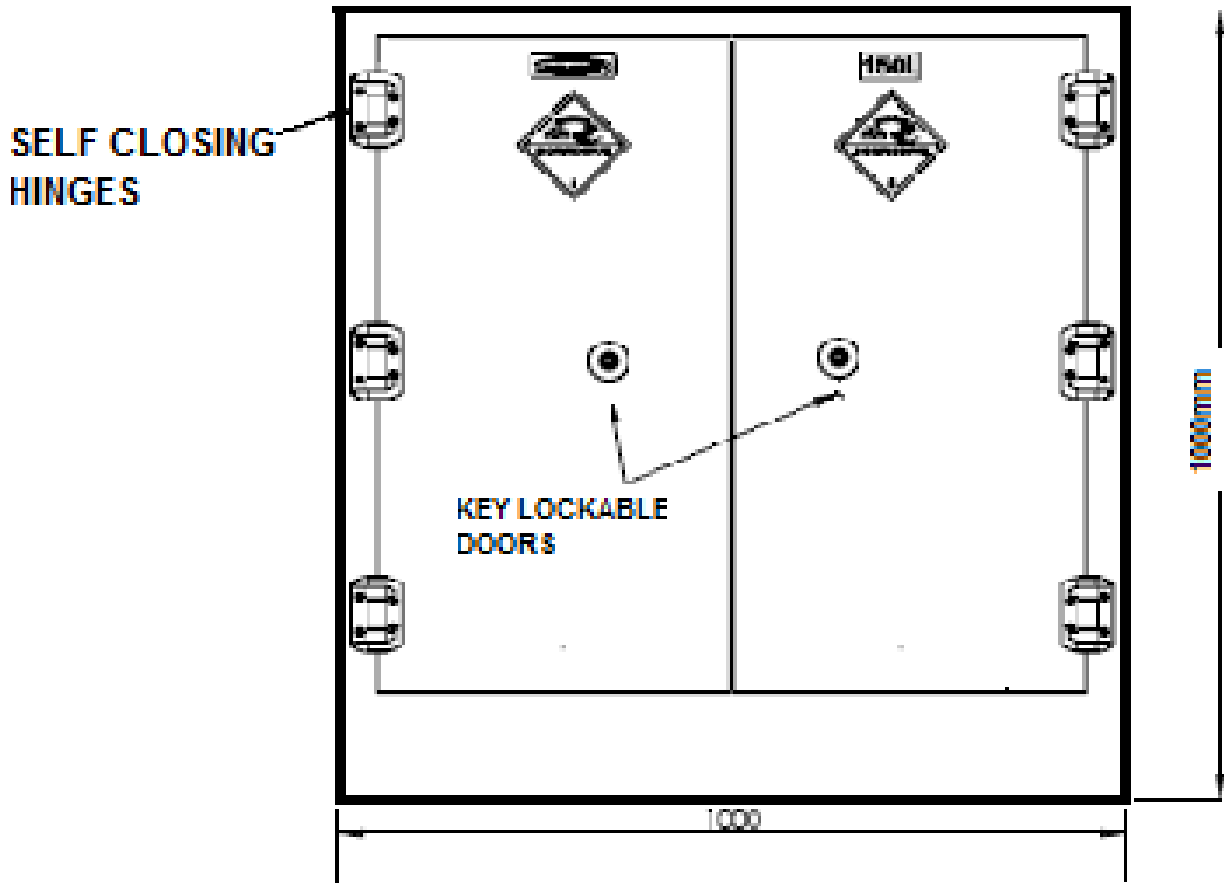

Side View

